



## Rhubarb


### Develosil ODS-5

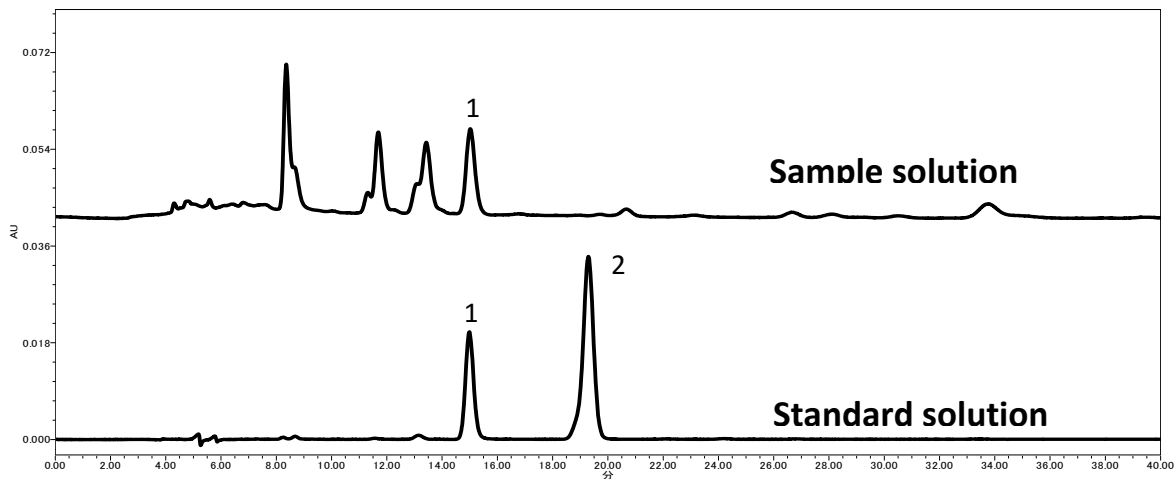

### Develosil ODS-UG-5

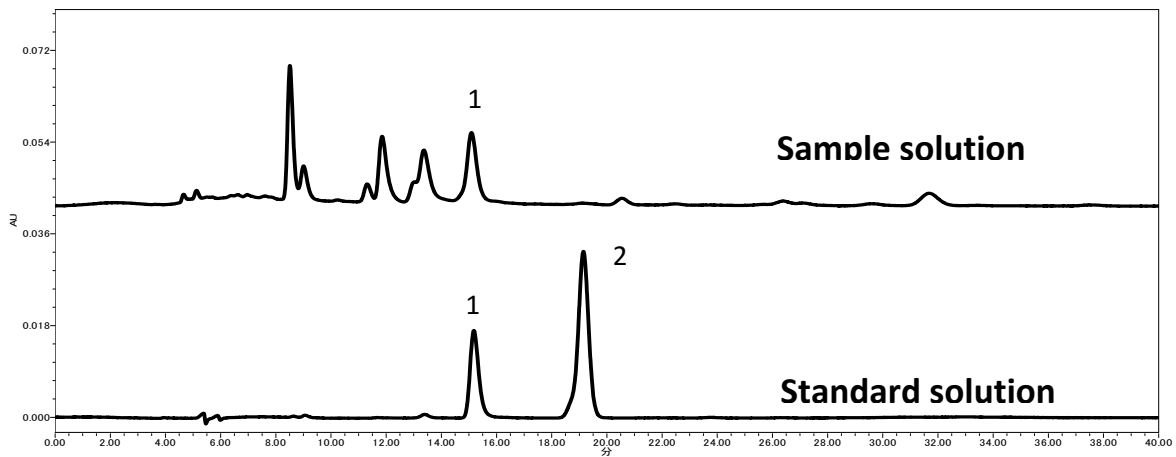

### Develosil XG-C18M-5

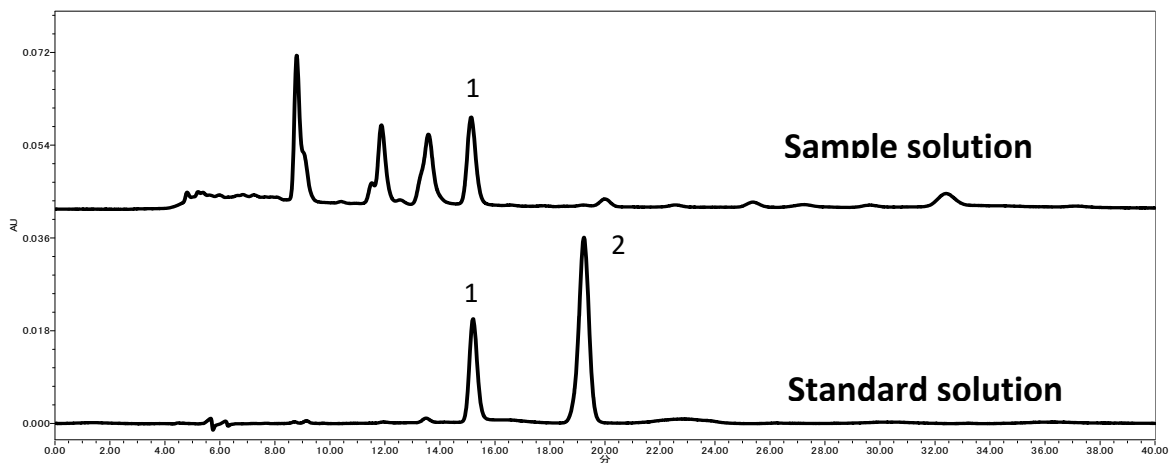

**Column: Develosil ODS-5 (4.6x150mm)**

**Develosil ODS-UG-5 (4.6x150mm)**

**Develosil XG-C18M-5 (4.6x150mm)**

**Mobile Phase: 1.25% Acetic acid/Acetonitrile=4/1**

**Flow rate: 0.34mL/min (Develosil ODS-5)**

**0.35mL/min (Develosil ODS-UG-5)**

**0.33mL/min (Develosil XG-C18M-5)**

**Temperature: 40°C**

**Sample: Standard solution**

**1: Sennoside A (0.05mg/mL)**

**2: Naringin (0.05mg/mL)**

**Injection Volume: 20µL**

**Detection: UV340nm**

	Requirement	Result
Resolution(1,2)	$\geq 3$	7.02 (Develosil ODS-5)
		6.36 (Develosil ODS-UG-5)
		6.69 (Develosil XG-C18M-5)