
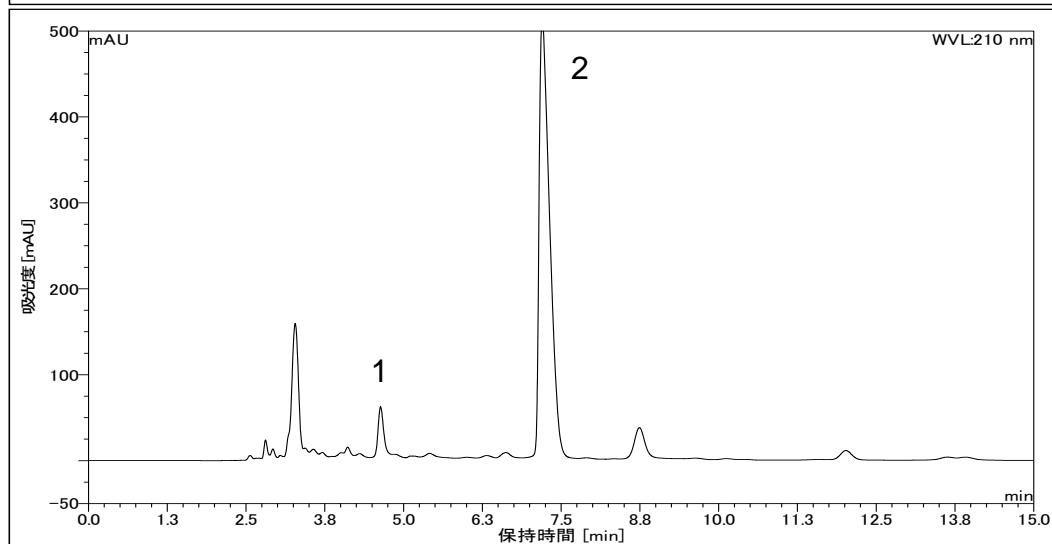


ドリンク中に含まれるアスコルビン酸・クエン酸の分析

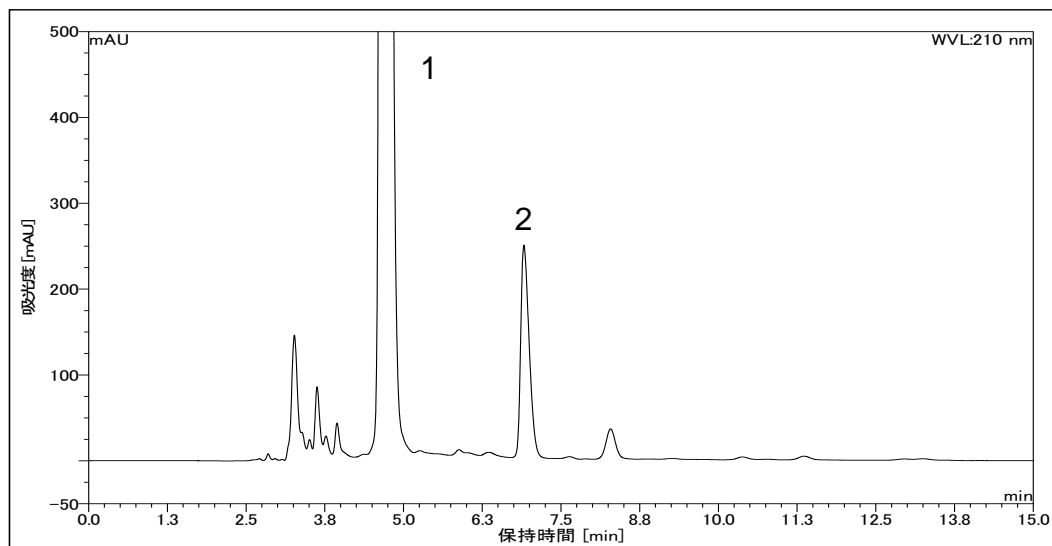

Conditions;

Column: Develosil XG-C1 (5um)

Size: 4.6x250mm

Mobile phase: Acetonitrile/0.2%Phosphoric acid=2/98

Flow rate: 1.0ml/min


Temperature: 30°C

Detection: UV210nm (PDA)


Sample: 1.Ascorbic acid 2.Citric acid

Injection volume: 10uL

System: DIONEX/Ultimate 3000 Series

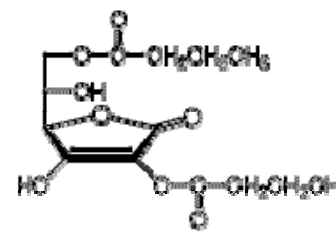

Ascorbic acid

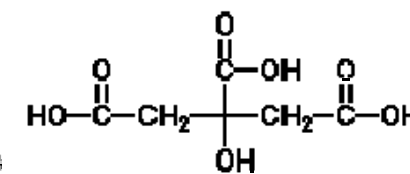

Citric acid