
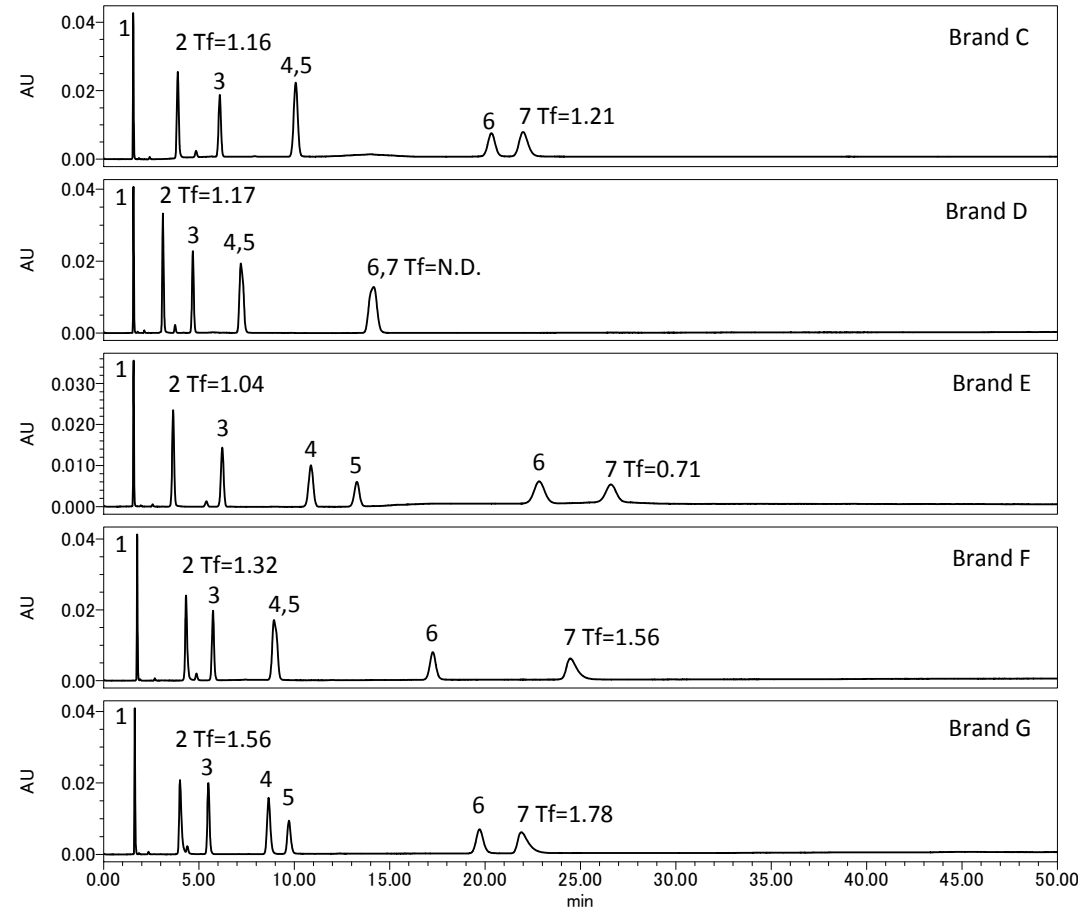


# Application Data

## No.092 Concurrent evaluation of ODS columns/ODSカラムの一斉評価

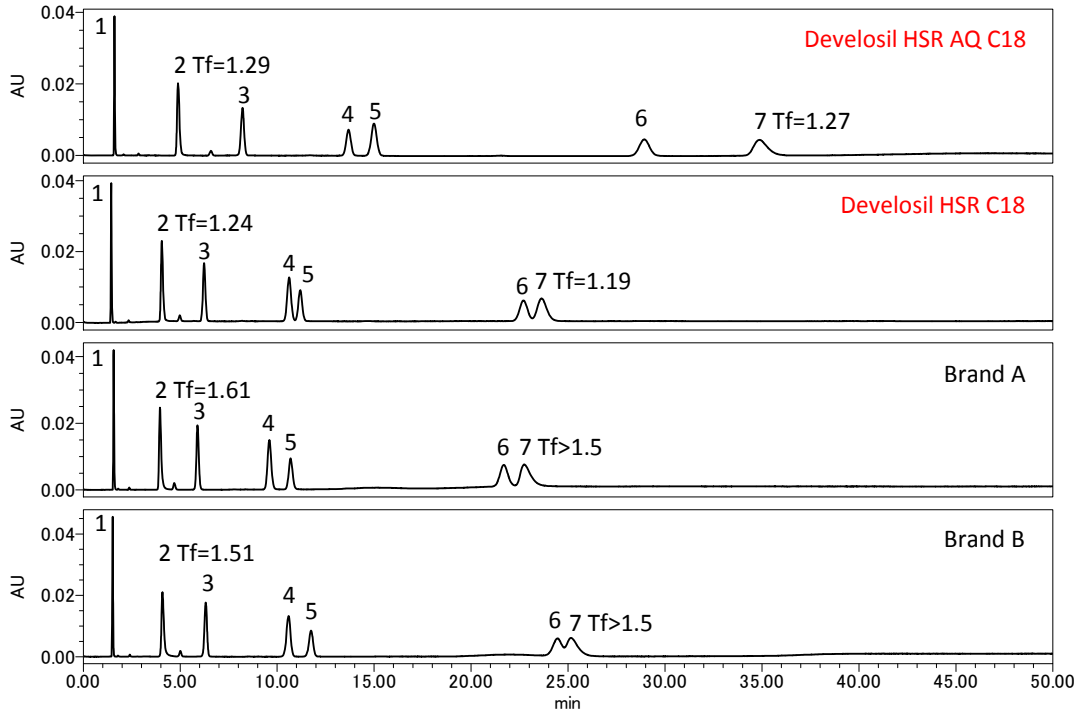

**Conditions;**

Column: C18 Column, 5um (4.6x150mm)  
 Mobile phase: Methanol/25mM Ammonium phosphate, pH7.0=65/35  
 Flow rate: 1.0mL/min Temperature: 40°C Detection: UV254nm  
 Sample: 1.Uracil 2.Propranolol 3.Butylparaben 4.Dipropylphthalate  
 5.Naphthalene 6.Acenaphthene 7.Amitriptyline  
 Injection volume: 2.0uL  
 System: Waters alliance

