


核酸塩基の分析

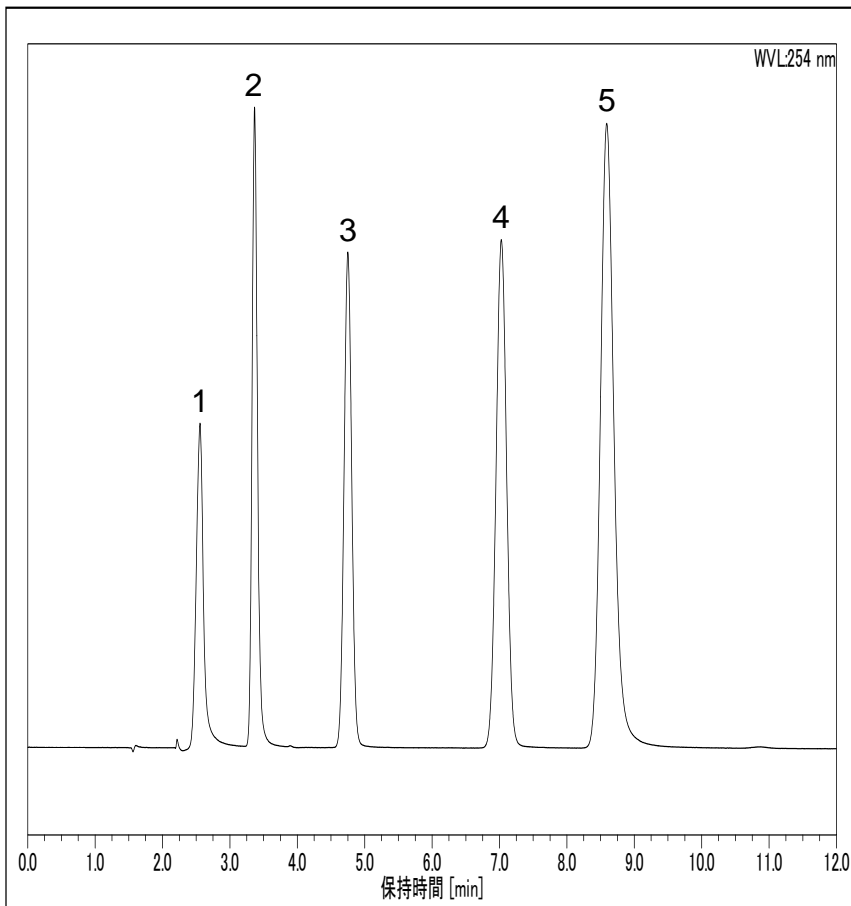

Sample

- 1.Cytosine (0.063mg/ml)
- 2.Cytidine (0.16mg/ml)
- 3.Uracil (0.062mg/ml)
- 4.Uridine (0.18mg/ml)
- 5.Thymine (0.20mg/ml)

Conditions;

Column: Develosil XG-C18M (5um) 4.6x150mm
Mobile phase: 50mM Phosphate Buffer (pH7.0)
Flow rate: 1.0ml/min
Temperature: 40°C
Detection: UV254nm
Injection volume: 2.0uL
System: DIONEX/Ultimate 3000 Series