


農薬の分析-3

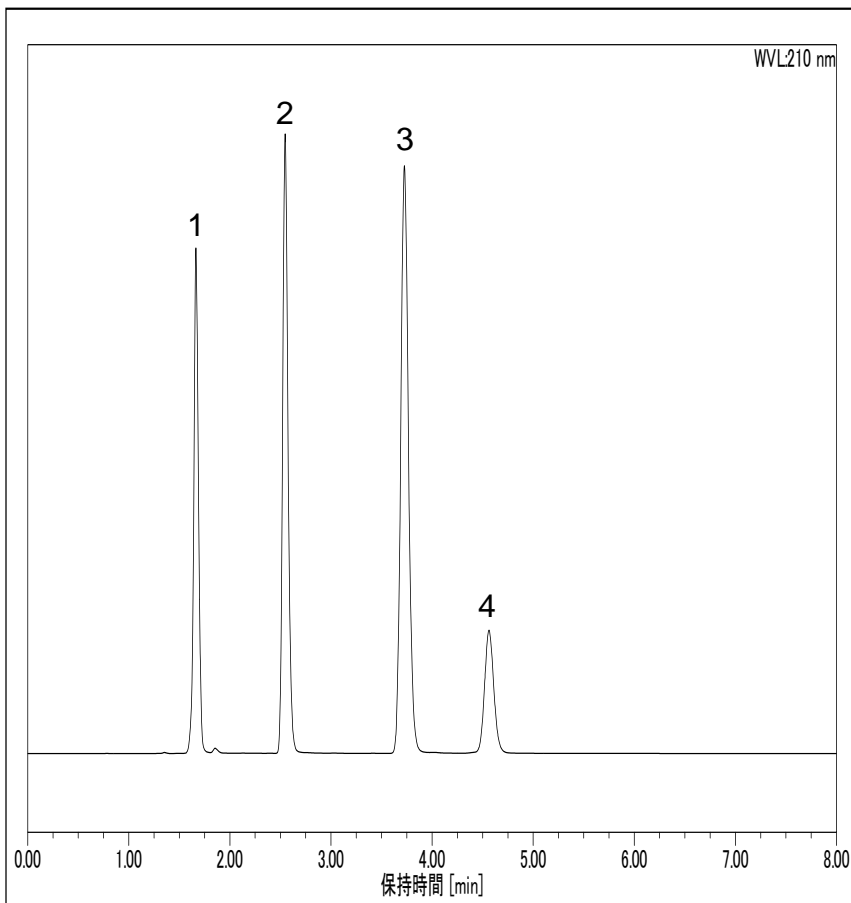

Sample

1.Asulam (0.41mg/ml)

2.Thiuram (0.36mg/ml)

3.Iprodione (0.39mg/ml)

4.Bensulide (0.42mg/ml)

Conditions;

Column: Develosil XG-C18M (5um) 4.6x150mm

Mobile phase: Acetonitrile/Water=70/30

Flow rate: 1.0ml/min

Temperature: 40°C

Detection: UV210nm

Injection volume: 2.0uL

System: DIONEX/Ultimate 3000 Series