


# Develosil Application Data

## No.059 Resveratrol in red wine/赤ワイン中のレスベラトロール

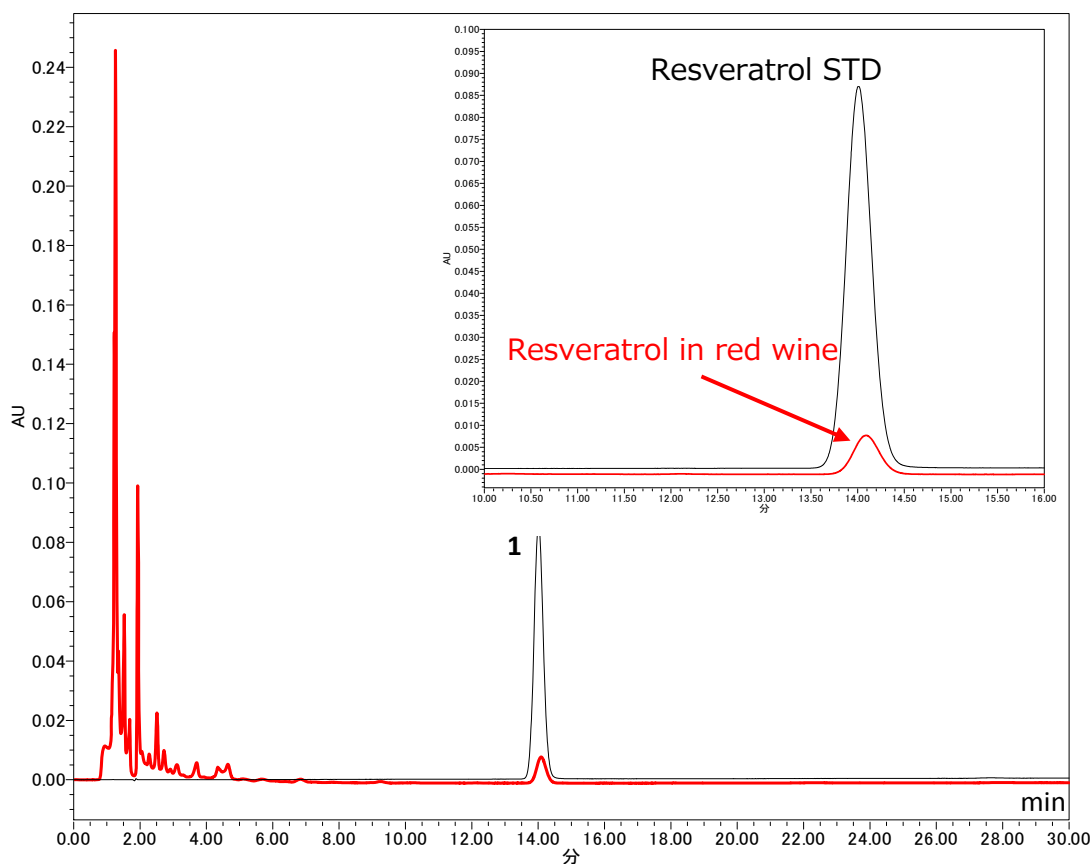

### Analytical Conditions:

Column: Develosil HSR C18,5 $\mu$ m (4.6x150mm)

Mobile phase: Acetonitrile/Water=20/80

Flow rate: 1.0mL/min

Temperature: 40°C

Detection: UV310nm

Sample: 1.Resveratrol (58ug/mL)

Injection volume: Resveratrol STD/5 $\mu$ L, Red wine/10 $\mu$ L

System: Waters alliance

※Red wine was filtered through a 0.45 $\mu$ m filter.

※赤ワインは0.45 $\mu$ mのフィルターでろ過しました。