


Develosil Application Data

No.083 Levofloxacin/レボフロキサシン

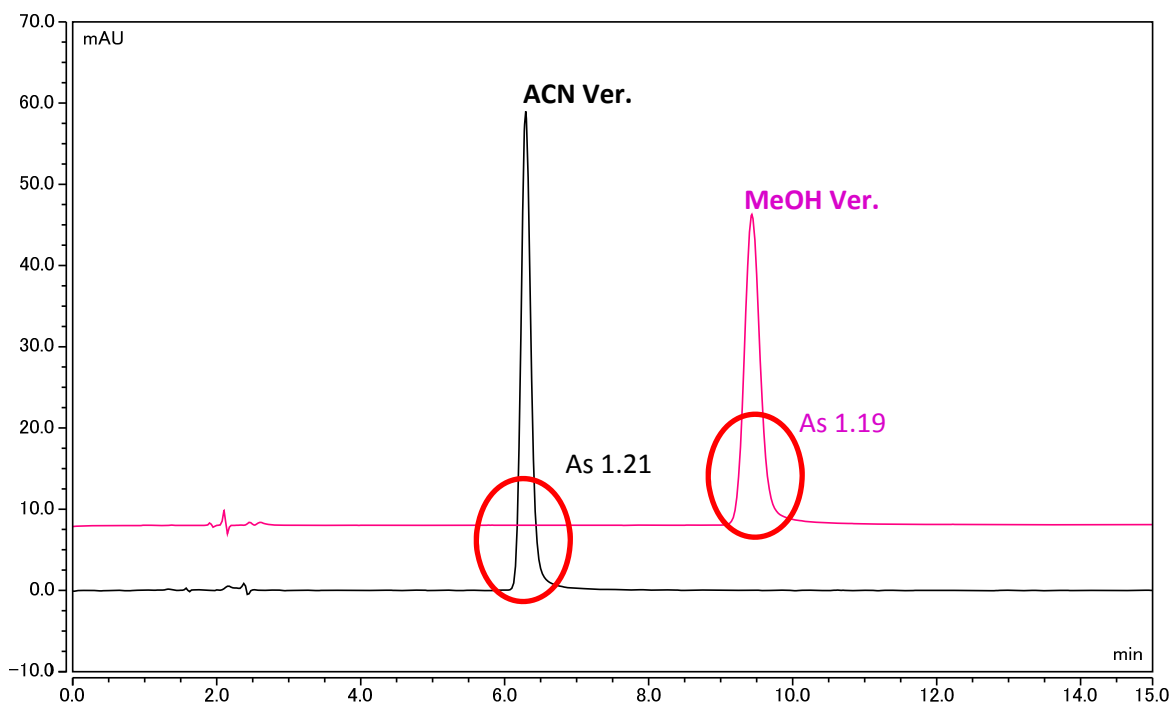

Analytical Conditions:

Column: Develosil HSR C18,5 μ m (4.6x150mm)

Mobile phase: Acetonitrile/25mM (NH₄)₂HPO₄, pH7.0 =20/80

Methanol/25mM (NH₄)₂HPO₄, pH7.0 =60/40

Flow rate: 1.0mL/min

Temperature: 40°C

Detection: UV294nm

Sample: 1. Levofloxacin (77ug/mL)

Injection volume: 2.0uL

System: Thermo Scientific Ultimate 3000