


Develosil Application Data

No.086 Peptide/ペプチド

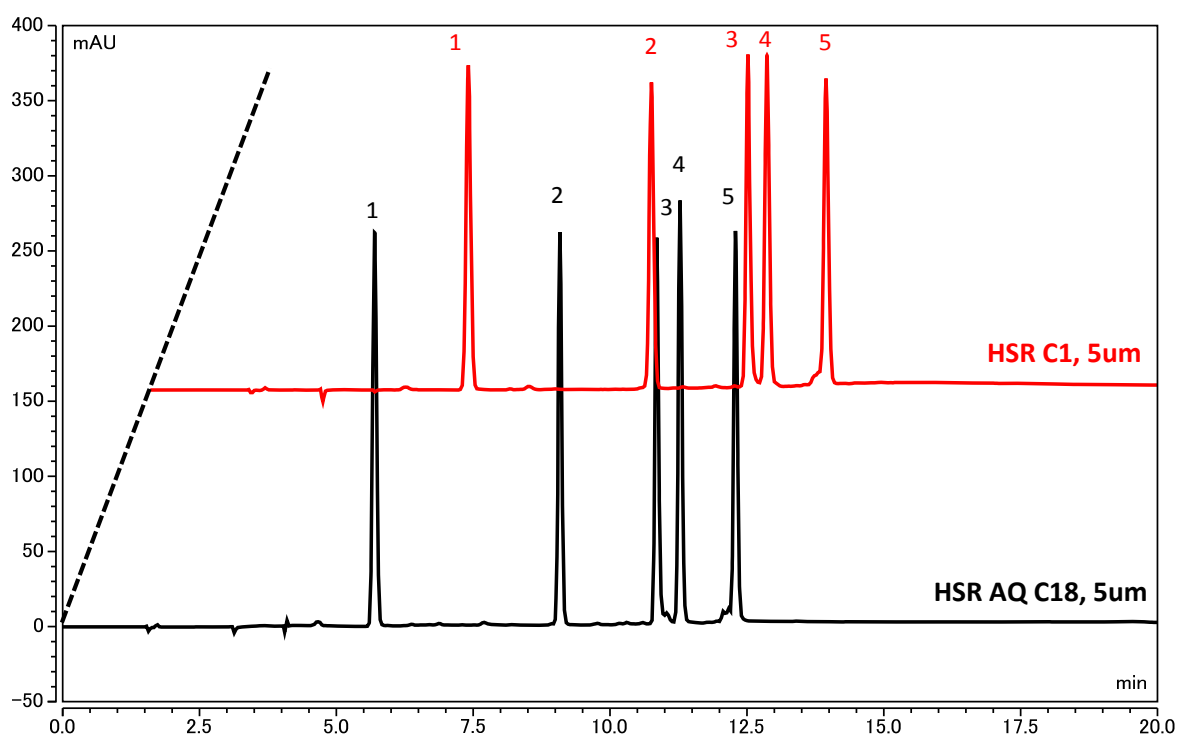

Conditions;

Column: **Develosil HSR C1, 5um (4.6x150mm)**

Develosil HSR AQ C18, 5um (4.6x150mm)

Mobile phase: A) 25mM NaCl (pH3.0) B) Acetonitrile

Flow rate: 1.0mL/min

Temperature: 30°C

Sample: 1. Gly-Tyr (0.05mg/mL)

2. Val-Tyr-Val (0.05mg/mL)

3. Met-Enkephalin (0.05mg/mL)

4. Angiotensin II (0.05mg/mL)

5. Leu-Enkephalin (0.05mg/mL)

Injection volume: 10.0uL

System: Thermo SCIENTIFIC Ultimate 3000

Gradient

Time (min)	Flow (mL/min)	%A	%B	Curve
0.0	1.0	100	0	5
10.0	1.0	70	30	5
20.0	1.0	70	30	5
20.1	1.0	100	0	5

